

NEXUS

Full-Arch Rehabilitation Redefined




The Complete Digital Workflow



Save Chair Time

Conventional protocols can take up to five patient visits. NEXUS iOS takes only two. Less chair time, less discomfort and accurate results.

NEXUS iOS offers a better alternative that results in an expertly engineered, industry-leading final prosthesis.

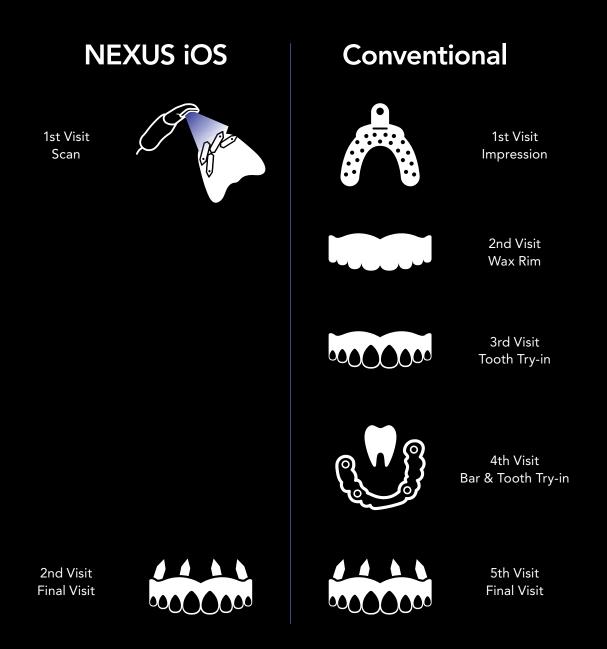

Up to 60% Reduction in Treatment Time

Potential of reducing three or more appointments from your treatment workflow.


What is NEXUS iOS?

A newly created first-of-it's-kind scanning and prosthesis production system.

Streamlining the development of custom, patient-specific full-mouth rehabilitation. Utilizing Osteon Medical's line of NEXUS prostheses.

Developed specifically for in-clinic use, from first scan to final insert.



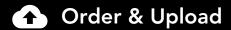







How it Works

With the use of NEXUS iOS Scan Gauge technology.

Scan Gauges allow for the capture of accurate intraoral scans that are used to create precise digital full-arch implant impressions.

No physical modeling needed. Compatible with all popular new and existing IO scanners, with a full suite of support offerings.



Scan Patient

Using your IO Scanner and NEXUS Scan Gauges.





Order & Upload

A new NEXUS prosthesis all through the NEXUS Ordering Platform.





Approve

The design of your case in 3D view with optional printed Try-In available.





Receive

Your finished NEXUS definitive case, ready for patient insert.



Lifelike Overlays

All NEXUS iOS cases are finished with expert technical artistry.

NEXUS iOS cases offer high-impact acrylic and monolithic multi-layered Zirconia finishing in a range of production methods to achieve true-to-life features like translucency. Finishing options and specifications can all be specified using the NEXUS ordering platform when uploading case data.

Our designers then provide a 3D file for digital approval, all contained on the cloud-based portal.



Lifetime Hardware Warranty

Using Osteon Medical's industry-leading implant-retained frameworks.

Osteon's titanium hardware underpins all NEXUS iOS prosthetics, with a lifetime warranty. Precision milling means precision patient fit every time.

Covering a range of indications, including:

- NEXUS Hybrid
- NEXUS Bridge
- NEXUS Removable

Limited Occlusal Space

The NEXUS Hybrid can be built to a minimum height of 6mm with monolithic zirconia, while maintaining significant strength.





NEXUS

Scan Gauges

The NEXUS Scan Gauge Kit and an IO scanner is all that is needed to begin performing scans for accurate full-arch rehabilitation.

The NEXUS Scan Gauge Kit is designed specifically for each clinician and their chosen abutment platform. Every kit has a unique identifying number and is associated with a registered NEXUS iOS account.

Each kit contains:

18 Scan Gauges in 3 lengths and 2 cuff heights

- 3 Scan Analogs
- 3 Narrow Scan Bodies

Driver

Unique Libraries

All NEXUS Scan Gauge Kits are designed to stringent engineering tolerances, produced by 5-axis milling machines. After production, exact dimensions of each Scan Gauge are measured and put into a custom library.

This is one of the unique steps that allows NEXUS iOS to ensure a high degree of accuracy.



Patient Journey

NEXUS iOS allows for a digital process throughout the entire patient journey.

is assigned to the clinician's

The NEXUS iOS workflow simply involves a sequence of scans and photos, uploaded to the NEXUS Ordering Platform at each stage of the patient journey.

Pre-Op

- Patient Photos
- Schedule Patient Surgery

Immediate Provisional

Surgery

account.

- IO Scans
- Download .STL
- Immediate Insert

Final Prosthesis

- Implants Integrated
- New IO Scans

Screw & Componentry Tray

- Order Final
- Final Prosthesis Insert



Global Headquarters

154 Middlesex Turnpike Burlington, MA 01803, USA Tel: +1 781-328-3490 Toll-free: 866-902-9272 www.KeystoneDental.com

Global KDG Offices:

USA

KeystoneDentalGroup Irvine 13645 Alton Pkwy Irvine, CA 92618, USA

Israel

Paltop Advanced Dental Solutions, Ltd. Hashita 5, Industrial Park P.O. Box 3568 Caesarea 3088900, Israel

Australia

Osteon Medical Headquarters 767 Springvale Road Mulgrave VIC 3170, Australia

France

Osteon Medical Europe 11 A, Bat B, rue des Aulnes Champagne au Mont d'or, 69410, France

Japan

Osteon Digital Japan 3-chōme-5-4 Nagayoshi Kawanabe Hirano Ward, Osaka, 547-0014, Japan

References

El-Haddad, H., Judge, R., Abduo, J. and Palamara, J., 2020. Laboratory Evaluation of Novel Implant Metal-Acrylic Prosthesis
Design: Influence of Monolithic Acrylic Veneer. The International Journal of Oral & Maxillofacial Implants, 35(1), pp.100-106.



